

Sterilisation Apparel & Accessories



High Performance Sterilisation Wrap
Innovation Focus: Barrier Performance



 **MULTIGATE**[®]

At Multigate we combine our deep understanding of how you work with our design capabilities to innovate and better serve you, your patients and our partners.

For 30 years we've created medical products, services and experiences that give you the advantage to perform at your best. We continue to do this.

Inspiring performances and better patient outcomes that will lead the world of hospital care to new standards.

We are Multigate.

We understand that sterile management, clinician protection and streamlined daily workflows are paramount in the CSSD.

With that in mind, we carefully designed our range of protective apparel, towels and tray liners and sterilisation accessories to power your performance.



MULTIGATE



Addressing CSSD clinician concerns about the time, costs and effectiveness of sterilising instruments, we introduced our simple to use, cost saving Pro-Zone® range of single use sterilisation wraps.

Made from polypropylene and using our 5 layer SSMMS protection technology, Pro-Zone® wraps offer an economic and time efficient sterilisation solution.



High Performance Sterilisation Wrap
Innovation Focus: Barrier Performance



TECHNOLOGY

Our Pro-Zone® range is made from polypropylene, utilizing state of the art 5 layer (Spunbond-Spunbond-Meltblown-Meltblown-Spunbond) technology.

The outer spunbond layers of the fabric provide strength, making it resistant to tears, punctures and abrasion, when correct storage and handling techniques are observed.


The inner layers of meltblown fibres provide a secure microbial and liquid barrier while still allowing effective sterilant penetration.

DESIGN

Our entire Pro-Zone® range of sterilisation wraps are a versatile high quality wrapping system that delivers many performance advantages including:

- Strength and resistance to punctures, tears and abrasions*
- Anti-static coating reduces static “shocks” to clinicians
- Low linting
- Soft and easy to fold
- Excellent drapability and low memory
- Excellent microbial barrier
- Effective sterilant penetration

* when correct handling techniques are used

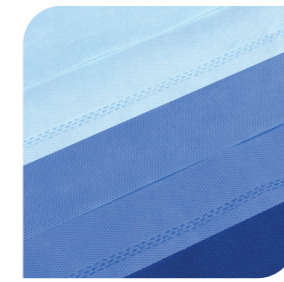

Pro-Zone®

The original Pro-Zone® sterilisation wrap.

Weights:

- 47gsm
- 61gsm
- 75 gsm

Compatible with:

- Steam
- ETO
- Low temperature gas plasma
- Low temperature steam/formaldehyde

Standards Compliance:

- AS/NZS 4187:2014
- ISO 11607



Pro-Zone®+2

Pro-Zone®+2 boasts the innovative use of 2 coloured sheets – blue and lilac, bonded together, enabling clinicians to rapidly identify any tears or damage in the wrap.

Weights:

- 47gsm
- 61gsm
- 75 gsm

Compatible with:

- Steam
- ETO
- Low temperature gas plasma
- Low temperature steam/formaldehyde

Standards Compliance:

- AS/NZS 4187:2014
- ISO 11607



Pro-Zone®+4G

Pro-Zone®+4G offers the protection of 2 bonded sheets, plus the added advantage of an absorbent layer to wick any residual condensate that may compromise the integrity of the wrap.

Weights:

- 61gsm
- 75 gsm

Compatible with:

- Steam
- ETO
- Low temperature steam/formaldehyde

Standards Compliance:

- AS/NZS 4187:2014
- ISO 11607

INFECTION PREVENTION

Pro-Zone®'s strong 5 layer construction reduces the risk of tears, punctures and abrasions, and delivers a low linting, high performing sterile barrier system.

- ✓ Infection prevention
- ✓ Process efficiencies
- ✓ Comfortable handling
- ✓ Space saving
- ✓ Recyclable
- ✓ Low linting
- ✓ Anti-static



Microbial barrier that delivers reliable barrier performance



Easy to fold material with low memory to provide excellent draping capability



PERFORMANCE

Streamlined Process.

Our Pro-Zone® range eliminates the costly decontamination processes associated with reusable sterilisation systems and delivers faster instrument processing.

The fabric is soft, has low memory and provides excellent draping capability allowing Pro-Zone® wrap to double as a sterile field, permitting the aseptic removal of pack contents and further streamlining theatre processes.

In-Service Support.

We understand that every CSSD is different. Our Sterilisation Specialists are trained in CSSD principles and practises. They are available to:

- Provide advice
- Deliver in-service training and support
- Evaluate sterilisation wrap performance.

FEATURES

Comfortable Handling.

Pro-Zone® wraps are anti-static, easy to fold, gentle on your hands and light offering a packaging system that is easy to handle.

Space Saving.

Pro-Zone®'s compatibility with a broad range of sterilants and its ability to cater for all instrument shapes and sizes with just a few different weights and specifications, reduces the number of sterilisation components and returns valuable space to CSSDs.

Recyclable.

The entire Pro-Zone® range is recyclable.

Pro-Zone®

The original Pro-Zone® sterilisation wrap.



| Pro-Zone® (47) Non-Sterile | | | COLOUR: BABY BLUE |
|----------------------------|-------------|--------------|---------------------|
| RE-ORDER | SIZE (mm) | WEIGHT (gsm) | CARTON QTY (Pieces) |
| 30-465 | 300 x 300 | 47 | 1000 |
| 30-464 | 450 x 450 | 47 | 500 |
| 30-460 | 610 x 610 | 47 | 500 |
| 30-470 | 710 x 710 | 47 | 300 |
| 30-461 | 915 x 915 | 47 | 150 |
| 30-467 | 1000 x 1060 | 47 | 150 |
| 30-466 | 1000 x 1400 | 47 | 100 |
| 30-462 | 1215 x 1215 | 47 | 100 |
| 30-468 | 1300 x 1525 | 47 | 100 |
| 30-463 | 1830 x 1370 | 47 | 100 |

| Pro-Zone® (61) Non-Sterile | | | COLOUR: BLUE |
|----------------------------|-------------|--------------|---------------------|
| RE-ORDER | SIZE (mm) | WEIGHT (gsm) | CARTON QTY (Pieces) |
| 30-471 | 710 x 710 | 61 | 200 |
| 30-472 | 915 x 915 | 61 | 150 |
| 30-475 | 1060 x 1200 | 61 | 50 |
| 30-476 | 1060 x 1400 | 61 | 50 |
| 30-473 | 1215 x 1215 | 61 | 50 |
| 30-479 | 1215 x 1215 | 61 | 50 |
| 30-477 | 1300 x 1525 | 61 | 50 |
| 30-478 | 1370 x 1370 | 61 | 50 |
| 30-474 | 1370 x 1830 | 61 | 50 |

| Pro-Zone® (75) Non-Sterile | | | COLOUR: DARK BLUE |
|----------------------------|-------------|--------------|---------------------|
| RE-ORDER | SIZE (mm) | WEIGHT (gsm) | CARTON QTY (Pieces) |
| 30-484 | 1060 x 1200 | 75 | 50 |
| 30-483 | 1060 x 1400 | 75 | 50 |
| 30-480 | 1215 x 1215 | 75 | 50 |
| 30-481 | 1370 x 1830 | 75 | 40 |

Pro-Zone®+2

Pro-Zone®+2 boasts the innovative use of 2 coloured sheets – blue and lilac, bonded together, enabling clinicians to rapidly identify any tears or damage in the wrap.



| Pro-Zone®+2 (47) Non-Sterile | | | COLOUR: BABY BLUE / LILAC |
|------------------------------|-------------|--------------|---------------------------|
| RE-ORDER | SIZE (mm) | WEIGHT (gsm) | CARTON QTY (Pieces) |
| 30-551Z | 300 x 300 | 47 | 480 |
| 30-552Z | 450 x 450 | 47 | 480 |
| 30-553Z | 600 x 600 | 47 | 240 |
| 30-554Z | 760 x 760 | 47 | 144 |
| 30-555Z | 910 x 910 | 47 | 144 |
| 30-556Z | 1010 x 1010 | 47 | 120 |
| 30-557Z | 1210 x 1210 | 47 | 48 |
| 30-558Z | 1370 x 1370 | 47 | 48 |
| 30-559Z | 1370 x 1820 | 47 | 48 |

| Pro-Zone®+2 (61) Non-Sterile | | | COLOUR: BLUE / LILAC |
|------------------------------|-------------|--------------|----------------------|
| RE-ORDER | SIZE (mm) | WEIGHT (gsm) | CARTON QTY (Pieces) |
| 30-560Z | 710 x 710 | 61 | 96 |
| 30-561Z | 910 x 910 | 61 | 72 |
| 30-562Z | 1010 x 1010 | 61 | 72 |
| 30-563Z | 1210 x 1210 | 61 | 24 |
| 30-564Z | 1370 x 1370 | 61 | 24 |
| 30-565Z | 1370 x 1820 | 61 | 24 |

| Pro-Zone®+2 (75) Non-Sterile | | | COLOUR: DARK BLUE / LILAC |
|------------------------------|-------------|--------------|---------------------------|
| RE-ORDER | SIZE (mm) | WEIGHT (gsm) | CARTON QTY (Pieces) |
| 30-570Z | 910 x 910 | 75 | 72 |
| 30-571Z | 1010 x 1010 | 75 | 48 |
| 30-572Z | 1210 x 1210 | 75 | 24 |
| 30-574Z | 1370 x 1370 | 75 | 24 |
| 30-573Z | 1370 x 1820 | 75 | 24 |



Pro-Zone®+2, 4G

Pro-Zone®+2, 4G offers the protection of 2 bonded sheets, the innovative use of 2 coloured sheets, blue and lilac, bonded together, enabling clinicians to rapidly identify any tears or damage in the wrap. plus the **added advantage of an absorbent layer** to wick any residual condensate.

| Pro-Zone®+2, 4G (61) Non-Sterile | | | |
|----------------------------------|--|--------------|---------------------|
| | | | COLOUR: BLUE |
| RE-ORDER | SIZE (mm) | WEIGHT (gsm) | CARTON QTY (Pieces) |
| 30-881 | 1010 x 1010 (absorbent layer 70 x 70) | 61 | 24 |
| 30-882 | 1210 x 1210 (absorbent layer 70 x 70) | 61 | 24 |
| 30-883 | 1370 x 1370 (absorbent layer 70 x 70) | 61 | 24 |

| Pro-Zone®+2, 4G (75) Non-Sterile | | | |
|----------------------------------|--|--------------|---------------------|
| | | | COLOUR: DARK BLUE |
| RE-ORDER | SIZE (mm) | WEIGHT (gsm) | CARTON QTY (Pieces) |
| 30-891 | 1010 x 1010 (absorbent layer 70 x 70) | 75 | 24 |
| 30-892 | 1210 x 1210 (absorbent layer 70 x 70) | 75 | 24 |
| 30-893 | 1370 x 1370 (absorbent layer 70 x 70) | 75 | 24 |
| 30-894 | 1370 x 1820 (absorbent layer 70 x 70) | 75 | 24 |



Towels and Tray Liners

Made from non-woven spun lace fabric cellulose/ polyester (MediClean™) & cellulose/polyester/ viscose (MediCEL®). These towels are ideal for a wide range of dry or wet CSSD applications. They are ideal for use as a tray liner and applications where the towel is soaked in water for long periods of time.

| MediClean Non-Sterile | | | | | |
|-----------------------|----------------------|---------|------------|--------------|---------------------|
| RE-ORDER | SIZE (cm) | COLOUR | PIECES/BOX | BOXES/CARTON | CARTON QTY (Pieces) |
| 57-112 | 31.8 x 33 (small) | ○ White | 100 | 12 | 1200 |
| 57-113 | 60 x 35 (large) | ○ White | 100 | 12 | 1200 |

| MediCel Non-Sterile | | | | | |
|---------------------|--------------------|---------|------------|--------------|---------------------|
| RE-ORDER | SIZE (cm) | COLOUR | PIECES/BOX | BOXES/CARTON | CARTON QTY (Pieces) |
| 57-120 | 35x30 (small) | ○ White | 100 | 12 | 1200 |
| 57-121 | 60 x 35 (large) | ○ White | 100 | 12 | 1200 |

Protective Apparel and Eyewear

- Thumb Secure Isolation/Protective Gown
- Liteback Thumb Secure Protective Gown
- Protective Eyewear

| Thumb Secure Isolation/Protective Gown Non-Sterile | | | | | |
|--|----------|---------|--------------|--------------|---------------------|
| RE-ORDER | COLOUR | SIZE | PIECES/BOXES | BOXES/CARTON | CARTON QTY (Pieces) |
| 29-285 | ● Yellow | Regular | 15 | 5 | 75 |
| 29-286 | ● Yellow | X-Large | 15 | 5 | 75 |
| 29-385 | ● Blue | Regular | 15 | 5 | 75 |
| 29-386 | ● Blue | X-Large | 15 | 5 | 75 |

| Liteback Thumb Secure Protective Gown Non-Sterile | | | | | |
|---|----------|----------|--------------|--------------|---------------------|
| RE-ORDER | COLOUR | SIZE | PIECES/BOXES | BOXES/CARTON | CARTON QTY (Pieces) |
| 29-291 | ● Yellow | Regular | 10 | 6 | 60 |
| 29-293 | ● Yellow | XX-Large | 10 | 6 | 60 |

| Protective Eyewear Non-Sterile | | | | |
|--------------------------------|---------|--------------|--------------|---------------------|
| RE-ORDER | COLOUR | PIECES/BOXES | BOXES/CARTON | CARTON QTY (Pieces) |
| 54-100NS | ● Black | 15 | 10 | 150 |



Thumb Secure Protective Gown



Liteback Thumb Secure Protective Gown



Protective Eyewear



Protective Corner

Our protective corner helps to reduce piercing risk while wrapping.

| Protective Corner Non-Sterile | | |
|-------------------------------|----------|---------------------|
| RE-ORDER | COLOUR | CARTON QTY (Pieces) |
| DP 2704 | ● Orange | 300 |

Empowering
performance
since 1986

Multigate Medical Products

27 Llewellyn Avenue

Villawood NSW 2163

Australia

T +61 (0)2 8718 2888

F +61 (0)2 8718 2889

enquiries@multigate.com.au

www.multigate.com.au

v13 - MAY 2024

For more information please contact your
local Business Development Manager:

(02) 8718 2888 or visit multigate.com.au